


LUTE +

Strontium based Glass Ionomer Luting Cement

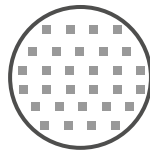

Ultra-Thin Film

Ultrafine particles (d50 – 5 microns) enable a thin film & Excellent Flow which is a vital requirement for a good luting cement



Extended Working Time




Cement remains in fluidic state up to 2 mins from start of mixing, giving the practitioner sufficient time for loading the cement and placement of Crown
**Please read IFU for mixing, working and setting time details*

Impermeable Glass Bottle

Powder is packed in a pharma grade Glass bottle with tamper evident cap

Accurate Dispensing

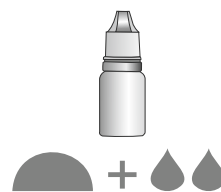


Dropper with metered dose enables accurate dispensing with ideal powder to Liquid ratio using 1 scoop of tapped powder & 2 drops of liquid



High Radiopacity

High Radio opacity allows for better Post-Operative Checks. This is due to the presence of Strontium in the Glass

Higher Radiopacity - Strontium Glass

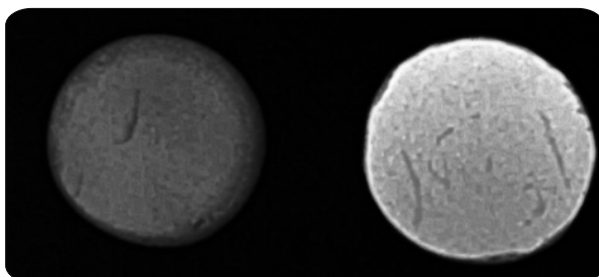

Packaging :

- Bottle of 12 g Lute Plus powder
- Bottle of 9 g Lute Plus liquid
- Mixing pad
- Spoon

(Shelf life: 2 years from date of manufacture)

