



# Universal Composite

NEW



*Best of both... handling & strength!*



High Radiopacity



Stackable Consistency  
Non Slumping

High Depth of Cure  
4-5 mm

Extremely high polishability  
& polish retention


Strength more than or equal  
to best in class composites



## Clinical Records

Case 1: Replacement of a non Anatomic Restoration  
- Dr. Srishti Punjani

Case 2: Bulk fill Restoration  
- Dr. Shushen Sharma

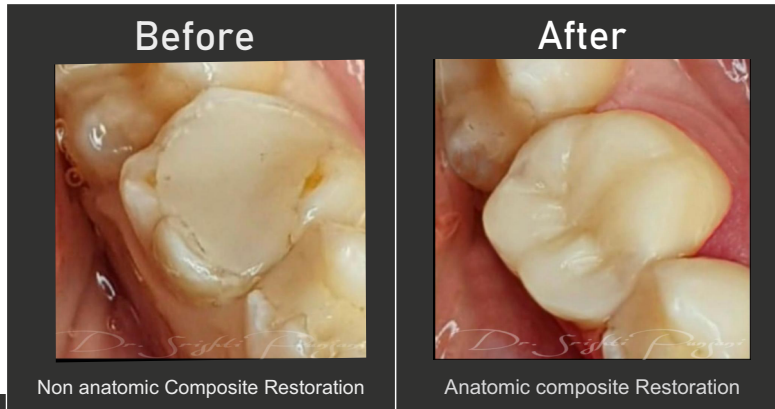

D Tech " Compoflo Universal " is a High Viscosity, Exceptional strength , Nano particle composite with shape holding , non - slumping properties. Some of the important features of " Universal" are :

- Our unique "PRIM technology" results in a depth of cure ranging from 4-5 mm with 20 seconds using a regular UV curing light allowing for its use as a bulk fill.
- Thanks to the well dispersed nano particles , it has exceptional polishability and polish -retention resulting in top notch, long term, aesthetic results.
- With Compressive strengths in the range of 300-350 Mpa & Flexural strengths in the range of 140-170 Mpa , this unique product exceeds or matches global leading highly filled composites.
- This product is a game changer and gives you best of both " **easy handling** like flowable and **high strength** like a high filled composite"... all combined into one easy to use product.




### PACKAGING

1 syringe of 2.5 gm universal composite  
2 nos dispensing tips

### Radiopacity comparison



Compoflo Universal

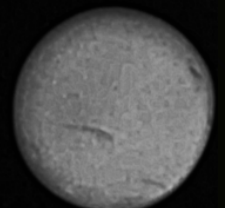

Leading Brand

**D- Tech Dental Technologies**  
62, Reality Warehousing  
Wagholi, Pune-412207 INDIA  
Url: [www.dtechasia.com](http://www.dtechasia.com)  
ig: @dtech3d

For orders contact :  
+91 9370145806  
+91 9370145807  
[info@dtechasia.com](mailto:info@dtechasia.com)

Proudly Made in India  
for the World



CE MDD Certified  
EN ISO13485:2016  
ICMED Issue 2  
REG NO.MFG/MD/2023/000135