
 D TECH DENTAL TECHNOLOGIES	Case Record		Case Date	06 Aug 2023
	Doctor Name : Dr.Priyam Aditya, Pune		Case Ref	Aug 23/DB
	Case No	CR-3DR-DB -02	Rev No &Date	NA

3D Printing - Partial Denture Case



Product used : Denture &C&B resin

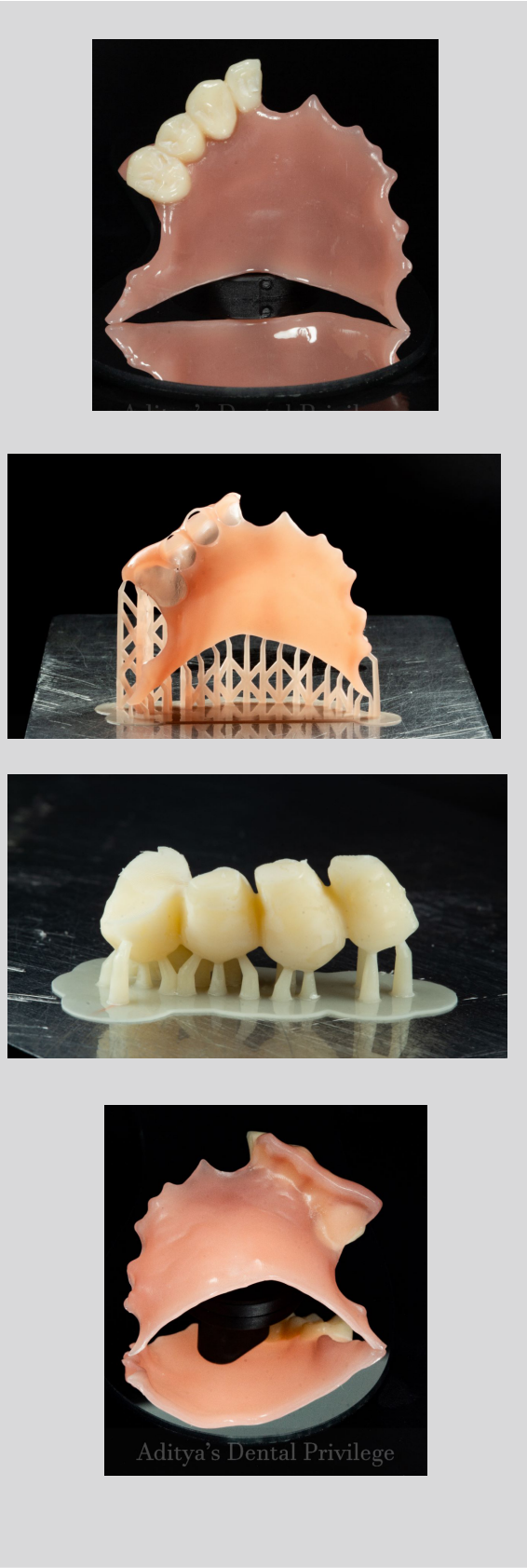

About The case
Case of 3D printing partial dentures

Outcome & Doctor's Comment
The patients oral mouth anatomy was digitally scanned and the stl file was imported into a denture designing software where one can begin to create a prosthetic. The Denture base was printed using Accuprint D tech Denture resin and the teeth were printed using C&B Accuprint resin. The two components were assembled and bonded together . The fitting was exceptionally good. The Denture material showed a high ability to be easily polished to a very highly aesthetic level.

About the Doctor
Dr. Priyam Aditya completed his B.D.S. and M.D.S. (Prosthodontics) from BharatiVidyapeeth University, Pune. He is member of Indian Dental Association (IDA), Indian Prosthodontic Society (IPS) and International Team for Implantolgy (ITI, Switzerland).

