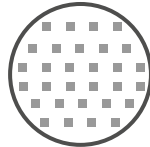


LUTE +

Strontium based Glass Ionomer Luting Cement

CE 0123

d-tech

Ultra-Thin Film

Ultrafine particles (d50 – 5 microns) enable a thin film & Excellent Flow which is a vital requirement for a good luting cement



Extended Working Time




Cement remains in fluidic state up to 2 mins from start of mixing, giving the practitioner sufficient time for loading the cement and placement of Crown

**Please read IFU for mixing, working and setting time details*



Impermeable Glass Bottle

Powder is packed in a pharma grade Glass bottle with tamper evident cap

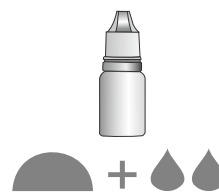

Accurate Dispensing


Dropper with metered dose enables accurate dispensing with ideal powder to Liquid ratio using 1 scoop of tapped powder & 2 drops of liquid



High Radiopacity

High Radio opacity allows for better Post-Operative Checks. This is due to the presence of Strontium in the Glass

Higher Radiopacity – Strontium Glass

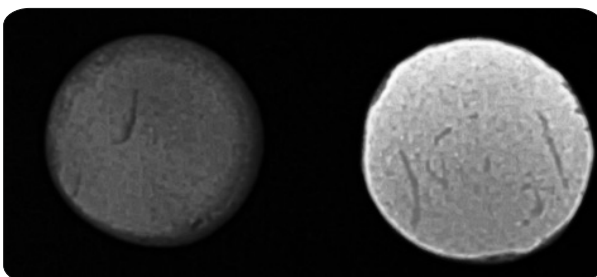

Packaging :

Bottle of 12 g Lute Plus powder

Bottle of 9 g Lute Plus liquid

Mixing pad

Spoon

(Shelf life: 2 years from date of manufacture)

D Tech Dental Technologies
Plot No 118/119, Sub Plot no 32
Lane No 7, Ramtekdi Industrial Estate
Hadapsar, Pune India 411028
Customer care- +91 9370145806
E mail: info@dtechasia.com
Url : www.dtechasia.com



CE MDD Certified
EN ISO13485:2016
ICMED Issue 2

Made in India
MFG LIC NO: MFG/MD/2024/000732



For orders contact :

+91 9370145806 +91 9370145807

+91 9370145808



info@dtechasia.com

ig: @dtech3d